

Anti GIP (Rat) Serum

Cat. No. Y102 Lot No. 655090129

Description: This antiserum was raised in a rabbit by immunization with a keyhole limpet hemocyanin (KLH) mixture of synthetic GIP (rat) peptide. The product vial contains 50 μ L of the titled antiserum obtained by lyophilizing its 0.001 M phosphate buffer (pH 7.0, 0.5mL) solution. It can be used for immunoassay, immunohistochemistry or other immunoreactions with GIP (rat).

Immunogen: Synthetic GIP (rat)/KLH mixture

Host: Rabbit

Amino Acid Sequence of GIP (rat)¹⁾:

YAEGTFISDY SIAMDKIRQQ DFNWLLAQK GKKNDWKHNL TQ

Product Form: Lyophilized unpurified serum

Size: 50 μ L

Reconstitution: Reconstitute the product with 0.5mL of 0.01M PBS (pH 7.0) to make a 10 fold diluted stock solution. If it is stored in a refrigerator, add moderate antiseptic to the solution (e.g. NaN₃ 0.1%).

Storage: The product will be stable for over one year if it be stored at -20°C to -80°C until opened. Upon reconstitution, the antiserum solution must be stored at 2°C to 8°C and used within one month. Repeated freezing- thawing should be avoided.

Suggested Working Dilution Range: 1:10,000~ for enzyme immunoassay (EIA); 1:5,000~32,000 for immunohistochemistry (frozen or paraffin sections). Optimal dilution should be determined by each laboratory for each application.

Specificity (based on non-competitive enzyme immunoassay): GIP (rat) 100%, GIP (18-42)(rat) 45%, secretin (rat) 0%, VIP (porcine) 0%, GLP-1 (7-36)-NH₂ <4.5%, GLP-1 (1-36)-NH₂ <3.8%, GLP-2 (rat) 0%, glucagon 0%.

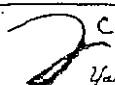
Positive Control (immunohistochemistry): Rat duodenum, jejunum, and ileum

Species Tested: Rat

REFERENCES:

- 1) Y. Higashimoto, J. Simchock, RA. Liddle, Molecular cloning of rat glucose-dependent insulinotropic peptide (GIP). *Biochim. Biophys. Acta.* 1132 (1):72-74, 1992

FOR RESEARCH LABORATORY USE ONLY
DO NOT USE ORGANIC SOLVENTS FOR DISSOLVING ANTISERUM


Yanaihara Institute Inc.

2480-1 Awakura, Fujinomiya, Shizuoka, 418-0011 Japan
FAX: +81-544-22-2770 TEL: +81-544-22-2771

For more information, please visit our internet website: <http://www.yanaihara.co.jp> or contact us to ask@yanaihara.co.jp

Copyright: Yanaihara Institute Inc. 2004 All rights Reserved